


Table Top Autoclave - Class B


Flash Autoclave

Robust and Compact Space-efficient

- B-process, N- process, Dri-Tec drying - process.
- Automatic (with manual override facility in case of auto mechanism failure).
- Sterilizing Cycle Processes available: 121 degree C Wrapped, 134 degrees C Wrapped, 134 degree C Flash/Rapid open instrument cycle, 134 degree C textile, Test programs: Bowie & Dick, Leak Test.
- Audio- visual Alarm facility available for notifying: Temperature sensor failure, Pressure sensor failure, Phase time out, Door improperly closed, Power failure (no alarm if duration of power failure < 10 sec), Safety device malfunctioning, Low water level, Cassette improperly closed.

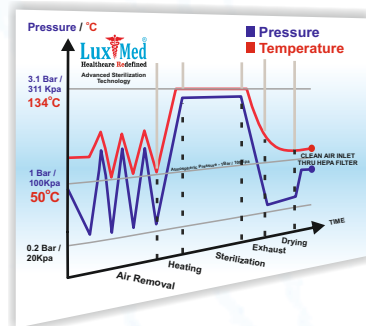

LuxMed
Healthcare Redefined

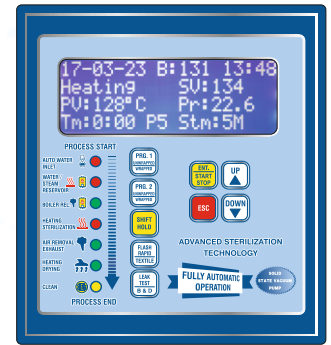


MODEL	CAPACITY
LM-TTA-15-K	15-18 LITER
LM-TTA-25-K	20-25 LITER
LM-TTA-25-G	20-25 LITER

Flash Autoclave Built In-Steam Storage

Fully Automatic Operation



- The tabletop autoclave can be used for medical and surgical items in ophthalmology, dentistry, and internal medicine clinics, such as packaged items, hollow and porous items for sterilization, and can also be used in emergency rooms and small laboratories.
- **LMA-ZB series sterilizer** is an automatic high temperature and rapid pressure sterilizer which works with steam as a medium.
- It can be widely used in departments of stomatology and ophthalmology, operating room, supply room, dialysis room, and other medical institutions. It is unpacked items, solid instruments, dental handpieces, endoscopes, implantable instruments, dressing fabric, rubber tubes, etc.
- **Build-in open-type water tank**
The sterilizer adopts an easy-clean open-type water tank which can support repeated program running if fully injected with water.
- **High-efficiency ultimate vacuum**
The sterilizer adopts a high-efficiency low noise vacuum system which has excellent effects.
- **Large LCD display for LMAZA/LMAZB**
The LCD screen can display temperature, pressure, time, operating status, failure warning, and other information.
It is convenient for customers to observe the sterilizer's running status.
- **Multiple program types**
The system has various programs, including packed items, unpacked items, BD testing program, vacuum testing program, and drying function. The custom program, rapid program and preheat function (for LMA-Z16B).
- **Standard USB port for LMAZA/LMAZB**
Users can store the sterilization data with a USB disk.
- **Optional mini printer can be attached to record the sterilization process.**

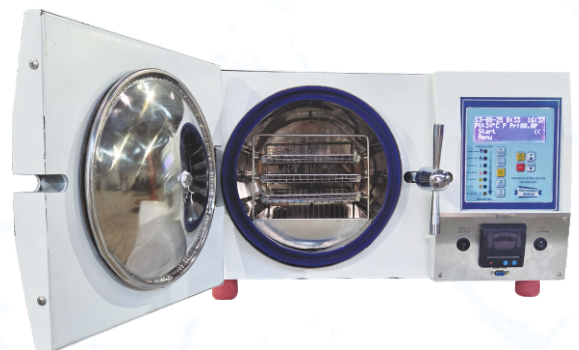

M: +91 8130383561, 9350831213 | E: cto@luxmed.in

SCAN & SEND MESSAGE :



SCAN & FOLLOW US ON :

